



Product List

IAS

Innovative Animal Supplies

Innovative Animal Supplies

Every Veterinary professional's trusted partner for transfusion solutions worldwide



IAS

Innovative Animal Supplies



TABLE OF CONTENTS


Canine Kits	Page 3-5	Mini 15	Page 10
Feline Kit	Page 6	CPDA1 Bag	Page 11
Large Animal Kit	Page 7	SA150 Blood Filter	Page 12
Transfer Kits	Page 8	Manual Press	Page 13
Veterinary Spike	Page 9	Contact us	Page 14

CANINE BLOOD COLLECTION KITS

Canine Blood Collection and Processing Kit (without leucodepletion filter)

Specific veterinary designed quad collection bag that allows collection of whole blood and processing into components including whole blood, packed red cells, plasma and platelet products.

Available in 300mL (>18kg donors) or 450mL (>25kg donors) collection volume options.

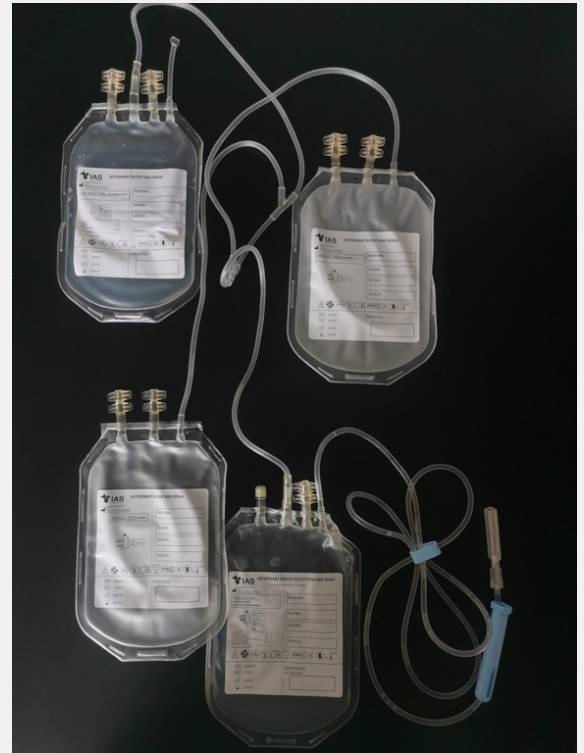

**Canine Quad Bag 300mL
Ref: DK4BWF2**


**Canine Quad Bag 450mL
Ref: DK4BWF1**

CANINE BLOOD COLLECTION KITS

Canine Blood Collection and Processing Kit (with leucodepletion filter)

Specific veterinary designed quad collection bag that allows collection of whole blood and processing into components including leukodepleted whole blood, leukodepleted packed red cells, leukodepleted plasma . Available in 450mL (>25kg donors).

Canine Quad Bag 450mL with filter
Ref: DK4BIF1

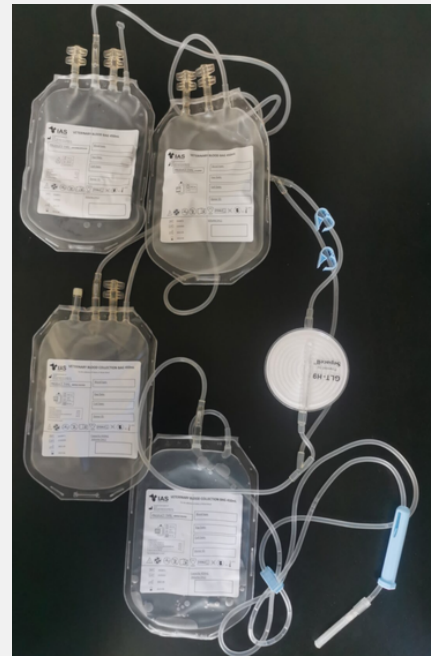

Canine Platelet Kit (with sampling pouch)

Specific veterinary designed quad collection bag that allows collection of whole blood and processing into components including whole blood, packed red cells, plasma and platelet products. Available in 450mL (>25kg donors).

Canine Platelet kit with sampling pouch
Ref: DK5BWF1



CANINE SINGLE COLLECTION KITS

Veterinary designed blood bag for the collection and storage of whole blood.

The bags are available in two sizes, 300mL (>18kg donors) or 450mL (>25kg donors).

May be used in canines or other veterinary species of a suitable weight and size.



Single Collection Kit 300mL
Ref: DK1BWF2

Single Collection Kit 450mL
Ref: DK1BWF1



FELINE BLOOD COLLECTION KIT

Closed system collection kit that includes CPD and SAGM and is designed for the collection of up to 50mL feline whole blood which may then be processed into red cells and plasma if desired.

The kit is supplied with a sterile 19 gauge needle with Luer for attachment at the point of use or may be substituted by the users own preferred needle.

The kit will support the collection of up to 100 mL of whole blood and may be useful for other small animal species.



Double Feline kit
Ref: CK2BWF1

LARGE ANIMAL COLLECTION KIT

Blood bag intended for the collection of up to 3L whole blood.

Anticoagulant of user's preference (e.g CPD, CPDA, ACD) can be added through any of the three connection ports and should be added at the manufacturers recommended volume for the calculated volume of blood to be collected from the donor.

12G siliconized sharp needle to ensure smooth venipuncture.

Long collection line allows tailored gravity flow rates onto a flat or suspended scale with the potential to allow rapid draw times for large donors and flexible location of suction units away from the donor.

No snap valves in collection lines reduces the chance of clots or sluggish flow within the line.



3L collection kit
Ref: CS1BWF1

TRANSFER BAGS

Designed for the transfer of up to 150 mL or 300 ml of blood components including packed red cells and plasma products.

May be attached to the parent pack using the port spike or alternatively the tubing is suitable for connection using most currently available sterile connection devices.



Transfer Kit 150mL
Ref: TK1BWF1

Transfer Kit 300mL
Ref: TK1BWF2

VETERINARY SPIKE

Unvented spike designed to fit all standard access ports on all collapsible blood and infusion bags. The spikes include a needleless access port with cover.




Veterinary spike
Ref: VS0BWF1

MINI 15 KIT

(Previously equine lysate kit)

Designed for the sterile separation of up to 15 aliquots of blood products and other fluids.

The Mini-15 kit contains 15 nested packs each capable of holding approximately 30 mL of fluid. The kit can be spiked on to your primary product pack or if preferred can be attached using a sterile connection device. For those wishing to limit the volume of product in each pack to 15 - 20mL a volume limiter device has been designed (see image right). IAS also provides dedicated storage boxes which can hold all 15 pouches either liquid or frozen. This kit has also been routinely used in the preparation and storage of equine platelet lysate.

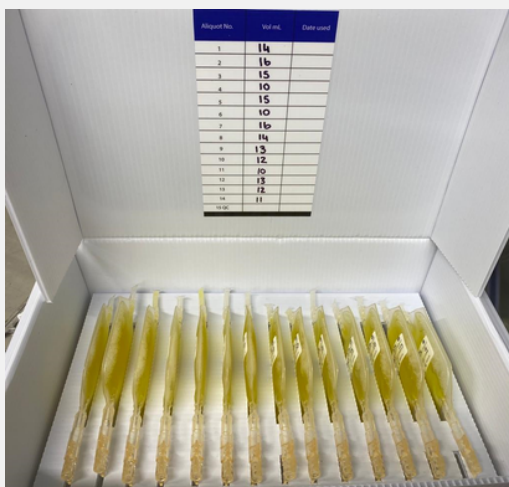

Equine Lysate kit
Ref: EL9BWF1

CPDA1 BAG


30mL of CPDA1 anticoagulant in a sterile pouch with both rubber needle access port and standard blood bag access port.



CPDA1 bag
Ref: CB1BWF1

150 MICRON AGGREGATE SYRINGE FILTER (SA150/BF150)

In-line filter for the removal of particulate material in blood components destined for transfusion into animals. Removes particles and aggregates of a size greater than 150 microns

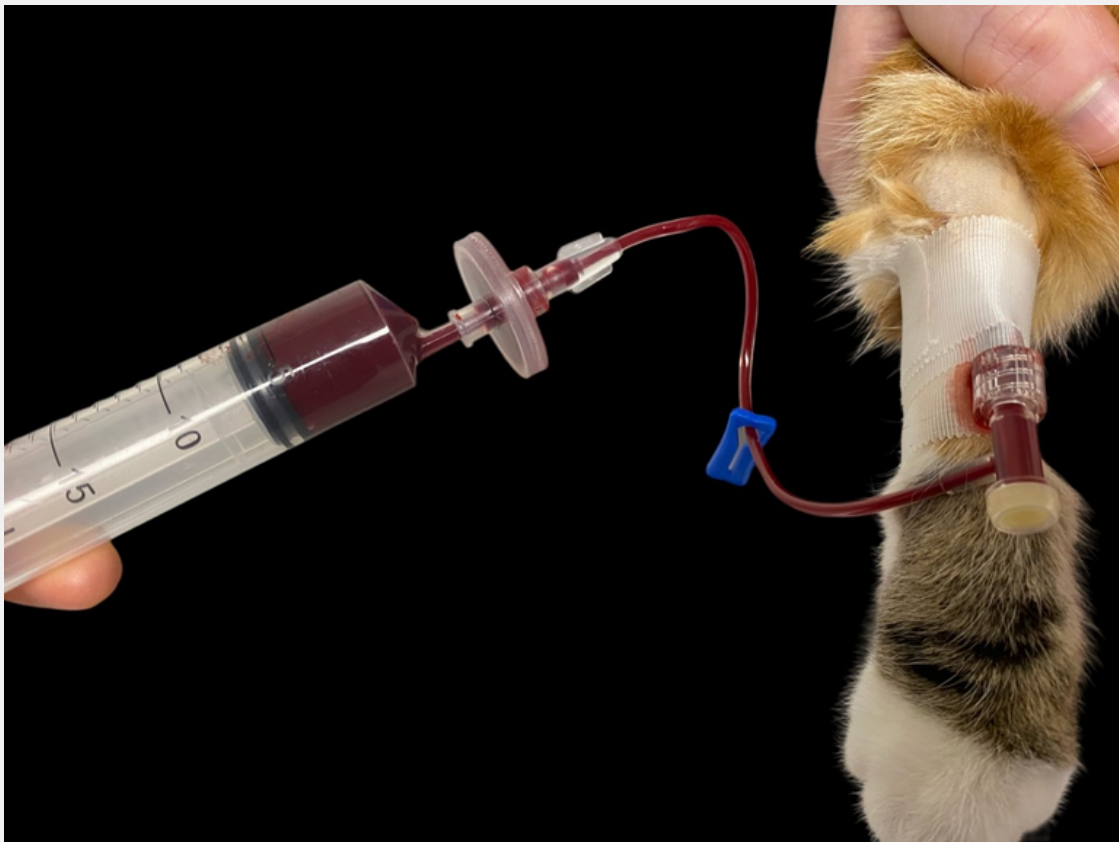

Syringe Filter
Ref: SA150

MANUAL PLASMA EXTRACTOR

The Manual Plasma Extractor is designed for the manual separation of blood components

Transparent, durable Plexiglass plate for easy visual monitoring of the process.

Lightweight and ergonomic design for user-friendly handling.

Robust construction ensuring long-term durability, even under intensive use.


Compatible with all standard blood bag types.

Made of corrosion-resistant materials, easy to clean and disinfect.

Provides uniform pressure across the blood bag for optimal plasma extraction.

Handle lock mechanism prevents unwanted movement during bag changes.

Ergonomic handle allows smooth and comfortable locking/unlocking of the device.

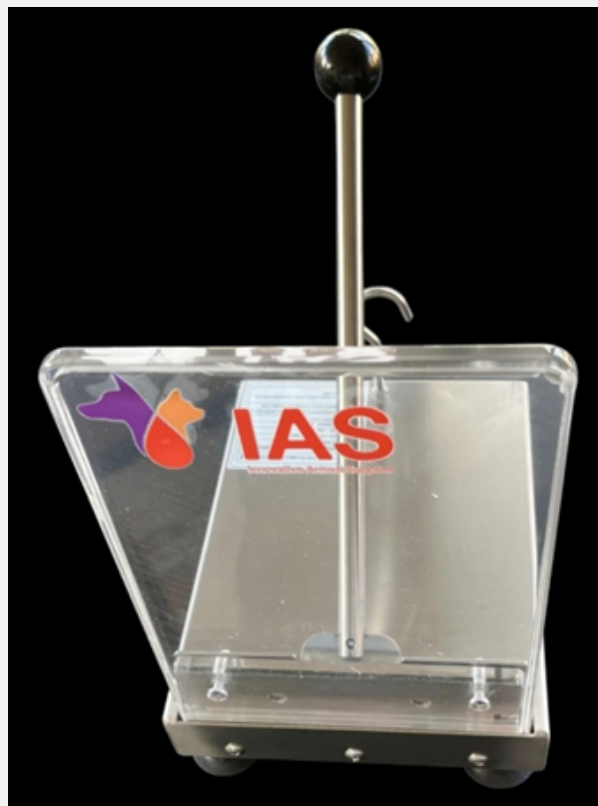

Ref: MP0BWF1

About us

We are five enthusiasts with a long and wide experience in the world of blood transfusion and medical consumable design and manufacture. Our backgrounds are in quality control, immunohaematology, veterinary blood banking and plastic engineering.

Our mission:

Every patient has the right to be transfused using the best available transfusion supplies which have been designed with the greatest consideration of individual species needs.

Our goal:

To assist veterinary professionals by providing safe and efficient blood collection and transfusion systems.

IAS Team



contact us

EMAIL

contact@innovativeanimalsupplies.com

WEBSITE

<https://innovativeanimalsupplies.com>

PHONE

Order Placement:
+44(0)7467917132

Product contact:
+44(0)7898890043

French contact:
+33(0)637897104

let's socialize!

Follow us on social media and don't forget to tag us in your photos!
[@innovativeanimalsupplies](#)